

Product Catalog

# Physiotherapy & Rehabilitation



**Compex FIXX 1.0 Massager**



# FIXX 1.0 MASSAGER

THE COMPEX FIXX 1.0 MASSAGER PROVIDES THE ULTIMATE DEEP TISSUE MASSAGE FOR SORE AND OVERWORKED MUSCLES

- Can help reduce muscle soreness
- Can help loosen tight muscles
- Can be used for warm up prior to activity

**2 YEAR WARRANTY**  
(EXCLUDING BATTERIES AND TIPS)



## SET CONTENT

- 1 CompeX Fixx 1.0 massager
- 1 premium travel case
- 1 battery • 1 power adapter
- 1 charging dock • 1 tip
- 1 instruction for use

## AVAILABLE ACCESSORIES

- 2 tips (cone and pad)
- Battery

<b>3 SPEEDS</b>	25sps, 33sps and 46sps*
<b>AMPLITUDE</b>	14mm
<b>WEIGHT (with battery and tip)</b>	990gr
<b>BATTERY LIFESPAN</b>	90min
<b>BATTERY CHARGING TIME</b>	120min

\*stroke per second



## HOW TO USE THE FIXX 1.0 MASSAGER


The CompeX® Fixx™ 1.0 massager provides the ultimate deep tissue massage for sore and overworked muscles. It can help reduce soreness, loosen tight muscles, and can be used for warm up prior to activity.



<b>TIPS</b>			
<b>USE</b>	<b>BALL</b>	<b>PADDLE</b>	<b>CONE</b>
	For general use	For large muscle group, IT bands	For small muscle and trigger points
<b>ARM POSITION</b>	0°, 45°, 90° position the arm as you like		
<b>SPEED</b>	<b>Speed 1: SOFT</b> For sensitive or sore muscles • <b>Speed 2: MEDIUM</b> For muscles that are not too small or sensitive <b>Speed 3: FIRM</b> For larger muscles that need something similar to a deep tissue massage		
	Keep moving slowly over muscles and avoid bony areas. Apply the pressure that is comfortable for you. Use on body parts for no more than 1 minute per specific muscle group. Combining speed and hand pressure customizes the massage to meet your needs and requirements for recovery.		

# HOW TO USE THE FIXX 1.0 MASSAGER

Refer to this chart to learn which muscle groups and body parts can benefit from the Fixx 1.0 massager.  
Try to follow the muscle fibers from origin to insertion and back again.  
Remember to use the device for both warm up and recovery to help take your training to the next level.

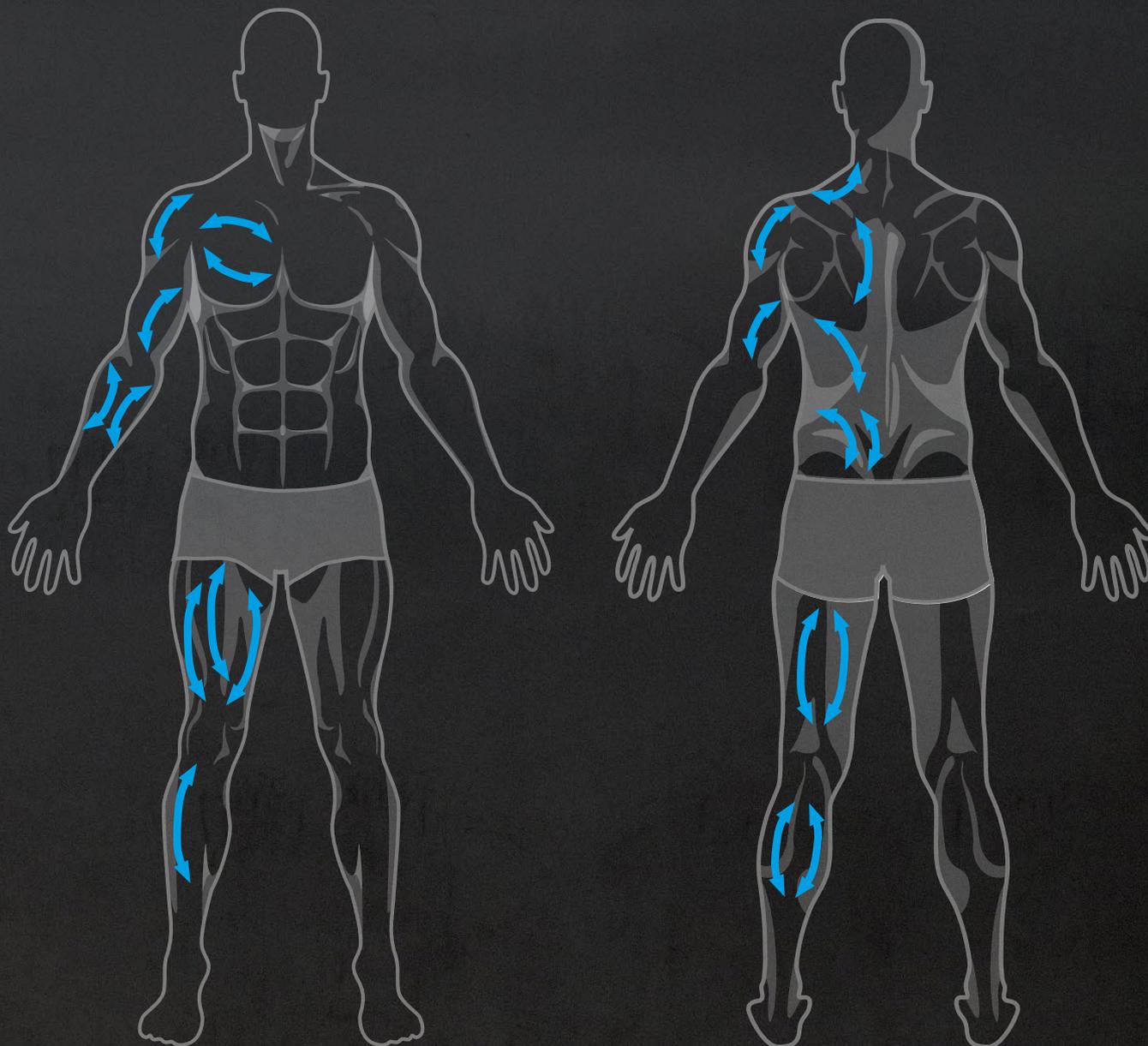

**COMPEX.COM**



**cnergy**<sup>®</sup>

**ONE STOP SHOP  
FOR ALL YOUR  
PHYSIO AND REHAB  
NEEDS**

Bracing & Supports  
Physiotherapy Modalities  
rPMS

Spinal Decompression Therapy  
Transcranial Magnetic Stimulation

Hot & Cold Therapy

Exercise Therapy

Cryo Chamber

Treatment Tables

Hydrotherapy Pool

Underwater Treadmill

Virtual Rehab

Balancing Platform

GAIT Lab

Isokinetic System

NeuroRehabilitation Robot

Robotic GAIT Trainer

Spine & Whole Body Analysis

Imported & Marketed by

**Cnergy Incorporation**

5277 | 8<sup>th</sup> Street | Ram Nagar North Extension | Madipakkam |  
Chennai 600091 | India, T: +91 44 4006 7974, +91 72000 53696 |  
incorporation@cnergy.in | customercare@cnergy.in | www.cnergy.in

Manufactured by  
**Enovis**  
USA